

KD-Validated Anti-Peroxiredoxin 1 Rabbit Monoclonal Antibody

Rabbit monoclonal antibody


ABG1065

Product Overview


Name	KD-Validated Anti-Peroxiredoxin 1 Rabbit Monoclonal Antibody
Accession(Primary)	Q06830
Host Species	Rabbit
Clonality	Monoclonal
Isotype	Rabbit IgG
Immunogen	A synthesized peptide derived from human PRDX1
Species Reactivity	Human, Mouse
Calculated MW	Predicted, 22 kDa , observed, 22 kDa
Storage Buffer	Supplied in PBS (pH 7.4) containing 50% glycerol, and 0.02% sodium azide.
Recommended Dilution	Western Blotting (WB): 1:1,000-1:5,000 Flow Cytometry (FCM): 1:2,000 Immunocytochemistry (IC): 1:1,000

Target information (Q06830)


Gene Name	PRDX1
About Gene Symbol	PRDX1 (Peroxiredoxin 1) is a peroxiredoxin antioxidant enzyme involved in oxidative stress response. It is an important research target in metabolic disease, cancer metabolism, and drug discovery and is widely studied to elucidate cellular mechanisms and identify therapeutic targets.

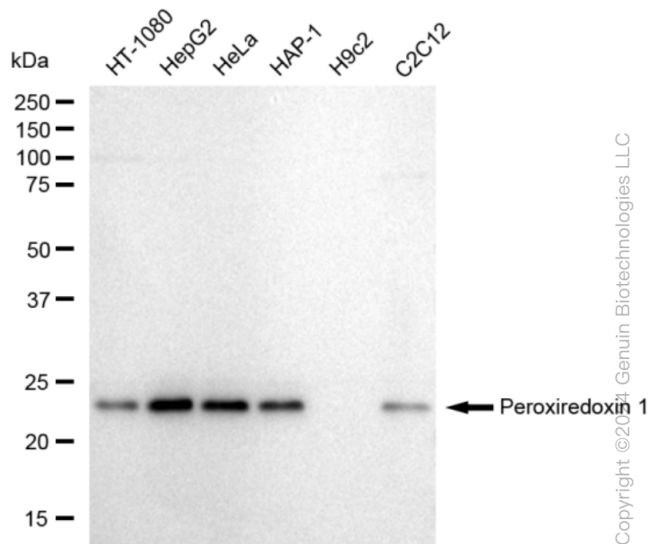


Western blotting analysis using anti-Peroxiredoxin 1 antibody (Cat#ABG1065). Total cell lysates (30 µg) from various cell lines were loaded and separated by SDS-PAGE. The blot was incubated with anti-Peroxiredoxin 1 antibody (Cat#ABG1065, 1:5,000) and HRP-conjugated goat anti-rabbit secondary antibody respectively.

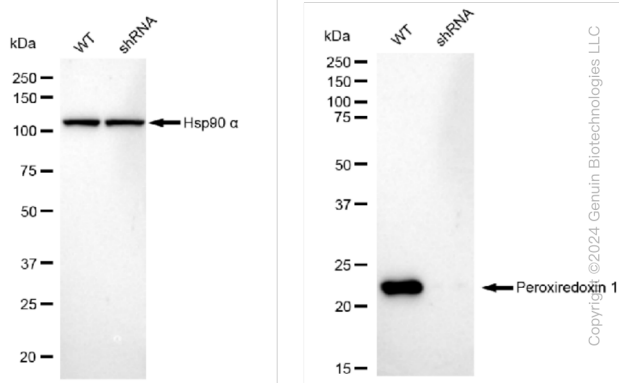

Western blotting analysis using anti-Peroxiredoxin 1 antibody (Cat#ABG1065). Peroxiredoxin 1 expression in wild type (WT) and peroxiredoxin 1 shRNA knockdown (KD) HeLa cells with 30 µg of total cell lysates. β-Tubulin serves as a loading control. The blot was incubated with anti-Peroxiredoxin 1 antibody (Cat#ABG1065, 1:5,000) and HRP-conjugated goat anti-rabbit secondary antibody respectively.

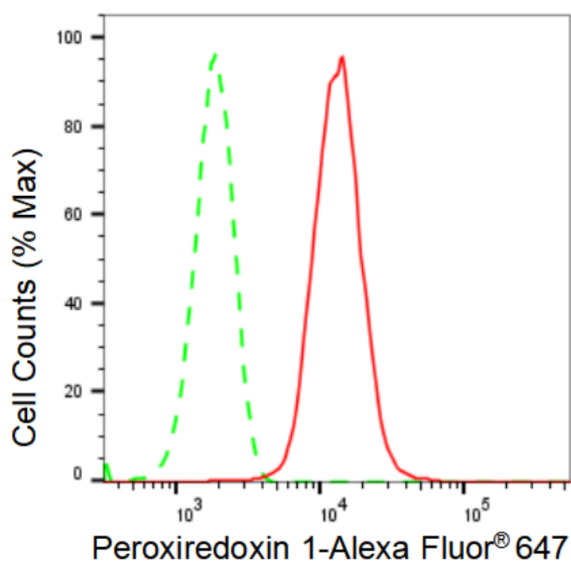

Flow cytometric analysis of Peroxiredoxin 1 expression in HepG2 cells using Peroxiredoxin 1 antibody (Cat#ABG1065, 1:2,000). Green, isotype control; red, Peroxiredoxin 1.

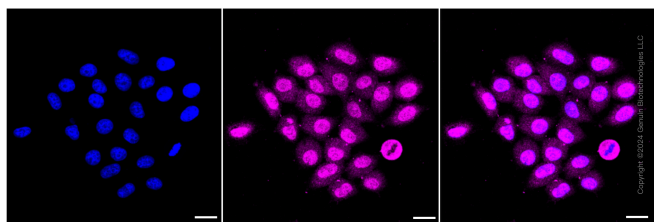

Immunocytochemical staining of HepG2 cells with Peroxiredoxin 1 antibody (Cat#ABG1065, 1:1,000). Nuclei were stained blue with DAPI; Peroxiredoxin 1 was stained magenta with Alexa Fluor® 647. Images were taken using Leica stellaris 5. Protein abundance based on laser Intensity and smart gain: High. Scale bar: 20 μ m.