

KD-Validated Anti-Nicastrin Rabbit Monoclonal Antibody

Rabbit monoclonal antibody


ABG1442

Product Overview


Name	KD-Validated Anti-Nicastrin Rabbit Monoclonal Antibody
Accession(Primary)	Q92542
Host Species	Rabbit
Clonality	Monoclonal
Isotype	Rabbit IgG
Immunogen	A synthesized peptide derived from human Nicastrin
Species Reactivity	Human, Mouse, Rat
Calculated MW	Predicted , 78 kDa , observed , 110-120 kDa
Storage Buffer	Supplied in PBS (pH 7.4) containing 50% glycerol, and 0.02% sodium azide.
Recommended Dilution	Western Blotting (WB): 1:1,000-1:5,000 Flow Cytometry (FCM): 1:2,000 Immunocytochemistry (IC): 1:1,000

Target information (Q92542)


Gene Name	NCSTN
About Gene Symbol	NCSTN (Nicastrin) is a protein kinase and signaling regulator involved in cell signaling and proliferation. It is an important research target in cancer, cell signaling, and therapeutic target research and is widely studied to elucidate cellular mechanisms and identify therapeutic targets.

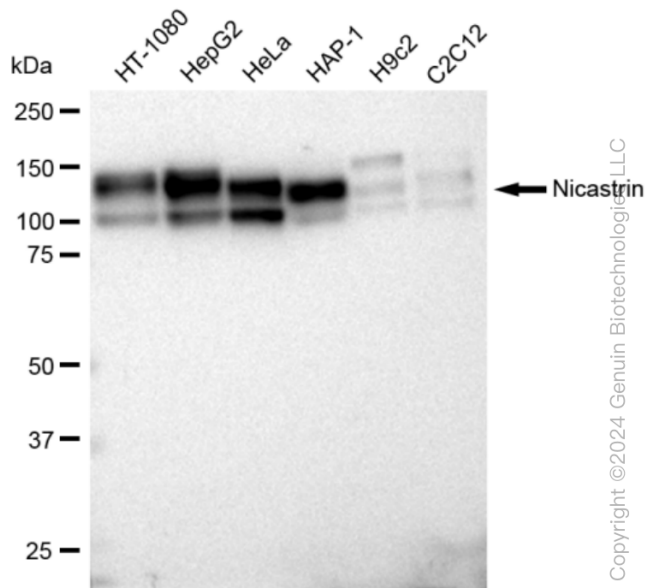


Western blotting analysis using anti-Nicastrin antibody (Cat#ABG1442). Total cell lysates (30 µg) from various cell lines were loaded and separated by SDS-PAGE. The blot was incubated with anti-Nicastrin antibody (Cat#ABG1442, 1:5,000) and HRP-conjugated goat anti-rabbit secondary antibody respectively.

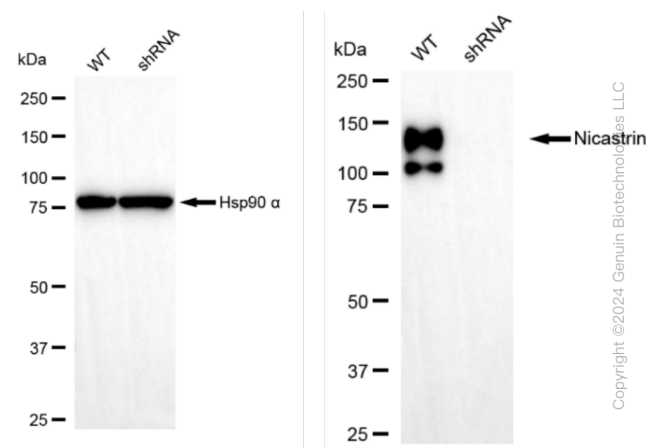

Western blotting analysis using anti-Nicastrin antibody (Cat#ABG1442). Nicastrin expression in wild type (WT) and Nicastrin shRNA knockdown (KD) HT-1080 cells with 30 µg of total cell lysates. β-Tubulin serves as a loading control. The blot was incubated with anti-Nicastrin antibody (Cat#ABG1442, 1:5,000) and HRP-conjugated goat anti-rabbit secondary antibody respectively.

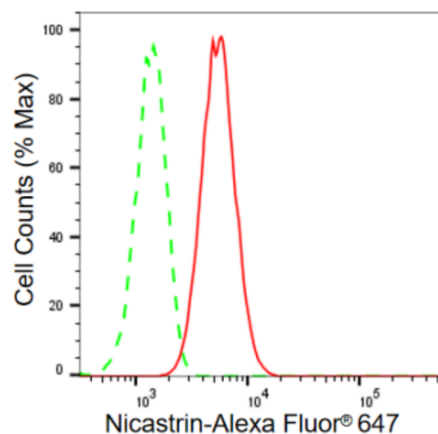

Flow cytometric analysis of Nicastrin expression in HepG2 cells using Nicastrin antibody (Cat#ABG1442, 1:2,000). Green, isotype control; red, Nicastrin.

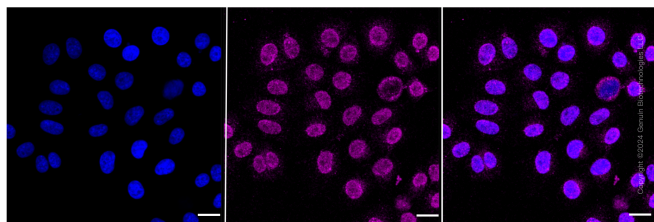

Immunocytochemical staining of HepG2 cells with Nicastrin antibody (Cat#ABG1442, 1:1,000). Nuclei were stained blue with DAPI; Nicastrin was stained magenta with Alexa Fluor® 647. Images were taken using Leica stellaris 5. Protein abundance based on laser Intensity and smart gain: Medium. Scale bar: 20 μ m.